
After installing Homer2 and AtlasViewer, test that the software is opening properly by following these steps. Make sure you don't get error messages, and that you see the screens displayed below.

Homer2


- Start Homer2:
 - If using the Homer2 executable: double-click the Homer2_UI icon on the Desktop (Homer2_UI.exe for Windows, Homer2_UI.command for Mac)
 - o If running Homer2 in Matlab: type Homer2 UI at the command line in Matlab

Note: *Homer2* can take several seconds to start.


• First, you will be asked to select a folder containing loadable .nirs files:

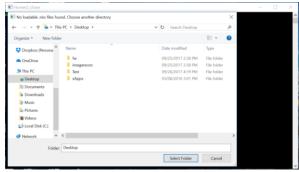

- If you're running the executable, navigate to Desktop > Test > Test-Homer2
- If running the Matlab version, navigate to PACKAGES > Test > Test-Homer2
 Select the Test-Homer2 folder:

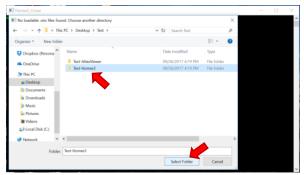


• Next you will be asked to choose a Process Options file (a .cfg file):

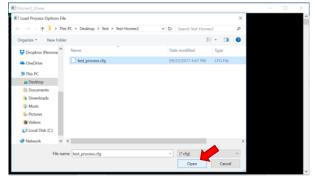


Select test_process.cfg and click Open


• The Homer2_UI window will open:

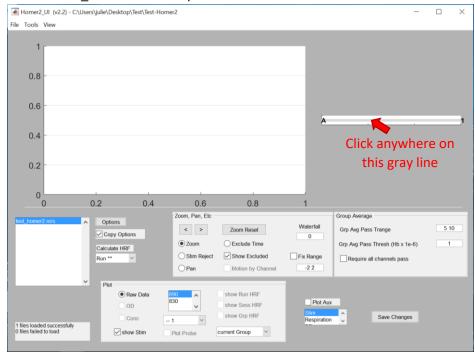

• Click anywhere on the grey line connecting (Source) A and (Detector) 1 on the right side of the window. You should see the following time series on the left:

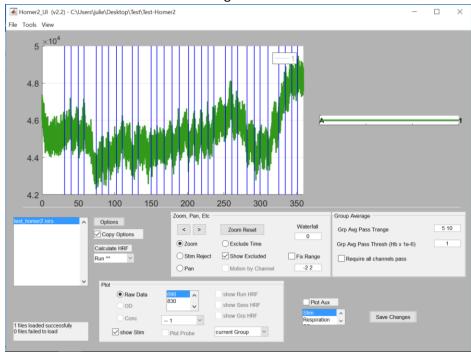


 Next select the 830 nm wavelength in the bottom window. You should see this new time series on the left:

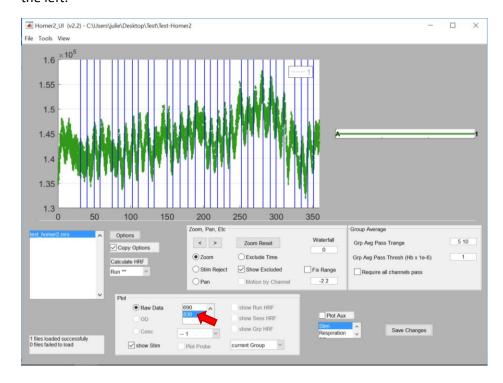


AtlasViewer

- Start AtlasViewer:
 - If using the AtlasViewer executable: double-click the AtlasViewerGUI icon on the Desktop (AtlasViewerGUI.exe for Windows, AtlasViewerGUI.command for Mac)
 - o If running *AtlasViewer* in Matlab: type AtlasViewerGUI at the command line in Matlab Note: *AtlasViewer* can take several seconds to start.
- First, you will be asked to select the subject folder:

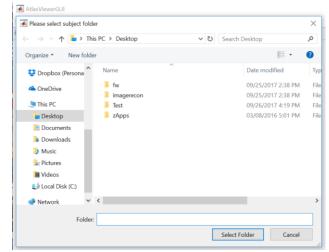

- If you're running the executable, navigate to Desktop > Test > Test-AtlasViewer
- If you're running Homer in Matlab, navigate to PACKAGES > Test > Test-AtlasViewer

Select the **Test-AtlasViewer** folder:

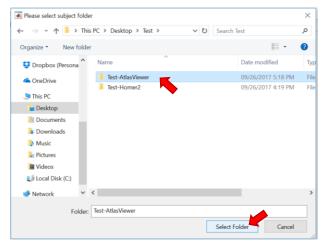

 The AtlasViewerGUI window will open and you should see a head with a probe on as displayed below:

